


ESD Results :

Size	Part No	Lot No	No. of trial	Decay Time (s)		Tribo Charge	SRM 100
				+ve 1 KV	-ve 1 KV		Ohms / sq
				<= 0.5 s		<= 20 Volt	<=10 ¹⁰
XS	TN1000W	N67 / HTPXS 8100006	1	0.1	0.2	6	3.55E+10
			2	0.1	0.1	10	3.91E+10
			3	0.1	0.1	8	2.57E+10
			Avg	0.1	0.1	8	3.34E+10
S	TN1001W	N67 / HTPS 8100021	1	0.2	0.2	17	2.54E+10
			2	0.2	0.1	12	2.02E+10
			3	0.2	0.1	15	3.19E+10
			Avg	0.2	0.1	15	2.58E+10
M	TN1002W	N67 / HTPM 8100023	1	0.1	0.1	9	2.55E+10
			2	0.1	0.1	6	1.96E+10
			3	0.1	0.2	10	2.74E+10
			Avg	0.1	0.1	8	2.42E+10
L	TN1003W	N67 / HTPL 8100021	1	0.1	0.1	14	2.52E+10
			2	0.1	0.2	16	2.99E+10
			3	0.1	0.1	12	2.02E+10
			Avg	0.1	0.1	14	2.51E+10
XL	TN1004W	N67 / HTPXL 8100014	1	0.2	0.1	14	3.77E+10
			2	0.2	0.2	8	3.92E+10
			3	0.1	0.2	12	3.04E+10
			Avg	0.2	0.2	11	3.58E+10
XXL	TN1005W	N67 / HTPXXL 8100006	1	0.2	0.2	9	3.47E+10
			2	0.2	0.2	13	3.66E+10
			3	0.1	0.2	11	3.01E+10
			Avg	0.2	0.2	11	3.38E+10

Conclusion

The Nitrile Glove White (Powder Free) under Batch No : N670710 passed the above specification.

Prepared by :

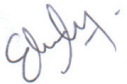

Ms Chong Ee Lane
QA Chemist

Approved by:



Ms Sabariah Salleh
QA Manager